



The Difference is **Precision.**



ASPIRATORS

Aspirators

Precision Medical Aspirators deliver top-of-the-line suction capabilities with outstanding safety features for a reliable medical device. Our range of portable and stationary aspirators are designed to carry out a wide range of suction procedures with continuous and intermittent operation. The Precision Medical line of aspirators feature hospital-grade power cords, replaceable inline hydrophobic bacteria filters, and hold multiple sizes of collection canisters, including 800cc, 1200cc and 2000cc.

EasyVac Aspirator (PM60)

Features

Lightweight at only 11.4 lbs

Durable design to withstand accidental knocks and bumps

Horizontal suction gauge for easy and comfortable reading

Suction cup feet securely hold the unit to a flat surface

Replaceable inline hydrophobic bacteria filter



Ordering Information

EasyVac Aspirator (PM60)


EasyVac includes 800cc canister, 6' patient tubing, 10" canister tubing, and replaceable inline hydrophobic bacteria filter

PM60 Accessories:



Precision Medical's EasyVac Aspirator provides reliable and safe suction of fluids, using up to 635 mmHg of vacuum. EasyVac provides additional safety to keep the suction canister from overflowing. The lightweight device is easy to transport and store, if necessary. The resetting thermal limit switch will protect the unit from overheating. The vacuum gauge is positioned horizontally for easy viewing from above. The gauge has a protective bracket to prevent unintended adjustments and knocks. Soft, pliable suction cup feet offer additional stability for the EasyVac Aspirator in any patient setting.

EasyGoVac Aspirator (PM65)

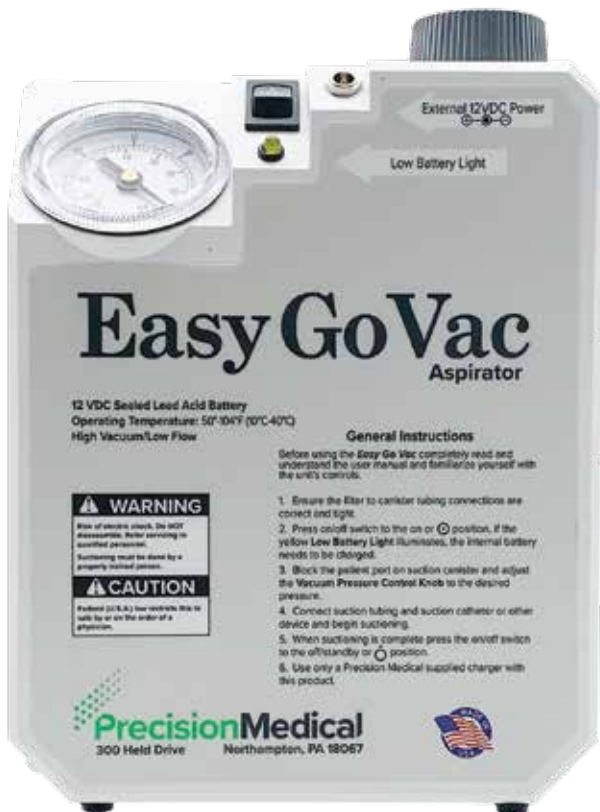

Features

Lightweight at only 9.6 lbs

2-hour battery runtime

Low battery light

Dual pump protection systems with float shut-off and replaceable inline hydrophobic bacteria filter safeguards

Can be used with any standard 800cc canister

Easy-to-read gauge

Ordering Information

EasyGoVac (PM65)

EasyGoVac includes 800cc canister, hospital-grade power supply and battery pack with charger, suction tubing, replaceable inline hydrophobic bacteria filter, elbow connector and nylon carrying case

PM65 Accessories:



The EasyGoVac battery-powered aspirator combines ease of use with powerful suction to create a portable unit that is ideal for hospitals, nursing homes, ambulances, and homecare. EasyGoVac produces vacuum up to 483 mmHg and provides dependable continuous suction power to remove respiratory fluids. The high-performance battery delivers 2 hours of runtime, more than any aspirator of comparable size. The low battery light lets users know when the unit needs to be recharged. The EasyGoVac is simple to clean and the closed design guards against bacterial growth. Safety features and safeguards in the pump prevent overflow and possible contamination.

PowerVac Aspirator (PM61)

ASPIRATORS

Features

Continuous operation

Automatic resetting thermal limit switch prevents overheating

Dual pump protection systems with float shut-off and inline hydrophobic bacteria filter safeguards

Easy-to-clean and infection control-friendly ABS plastic cabinet

Replaceable inline hydrophobic bacteria filter

Ordering Information

PowerVac Aspirator (PM61)

Includes PowerVac with 2000cc canister, 6' patient suction tube, 10" suction canister tubing, inline hydrophobic bacteria filter and hospital-grade power cord

PM61 Accessories:



The Precision Medical PowerVac Aspirator is built around a high-performance suction pump that delivers up to 622 mmHg of vacuum. The wide range of suction power allows for operation in homecare and hospital settings. PowerVac is designed for continuous use and has a built-in thermal safety system that protects the device from damage due to overheating. It features a dual protection system that includes a float shut-off in the canister lid and an inline hydrophobic bacteria filter to prevent back flow of fluids and some airborne contaminants from entering the pump. The PowerVac also has a unique, movable canister holder which allows the caregiver to rotate the canister to either side of the unit for better placement near the point of care.

PowerVac+ Aspirator (PM63)



Features

Continuous and intermittent operation options

Hospital grade vacuum regulator

Self-lubricated oil-less bearings

Movable canister holder for placement on either side

Collection canister with overflow protection

Easy-to-clean and infection control-friendly ABS plastic cabinet

Replaceable inline hydrophobic bacteria filter

Ordering Information

PowerVac+ Aspirator (PM63)

PowerVac+ includes 2000cc canister, 6' patient suction tubing, 10" suction canister tubing, replaceable inline hydrophobic bacteria filter, plastic elbow connector and hospital-grade power cord

PM63 Accessories:



The Precision Medical PowerVac+ Aspirator is designed for continuous and intermittent use. This powerful and durable suction pump combines the performance expected in a hospital with the ease of operation needed for homecare use. It features a 0-200 mmHg adjustable vacuum regulator gauge and delivers maximum vacuum of up to 250 mmHg. The PowerVac+ is a simple-to-use stationary aspirator with a durable design to withstand the rigors of hospital and homecare use. It features a hospital grade vacuum regulator, an inline hydrophobic bacteria filter to prevent contamination, self-lubricated oil-less bearings, and includes a collection canister with overflow protection.

Ways to place an order:



Contact your Precision Medical sales representative



Visit precisionmedical.com



Contact Customer Service:
orders@precisionmedical.com
1-800-272-7285

Do you need a special configuration?

Custom configurations, international specifications, special ratings and color options are available.

Contact Customer Service at 1-800-272-7285 for custom requests and more information.